


Hydrophobic coating
allowing for effortless & durable gliding
properties, no matter if wet or dry.


🎯 Rethink the basics
Passionate about bringing innovative ideas and medical device solutions to life, we are a team building on nearly 20 years of experience in the design and manufacturing of medical device solutions for partners around the world.
Even small changes in product design and quality can make a clear difference – making procedures easier and simpler with better outcomes.
Continually striving for solutions to the most basic problems, we don't go back to basics, we rethink the basics.



Atraumatic tip
designed for smooth and easy insertion.

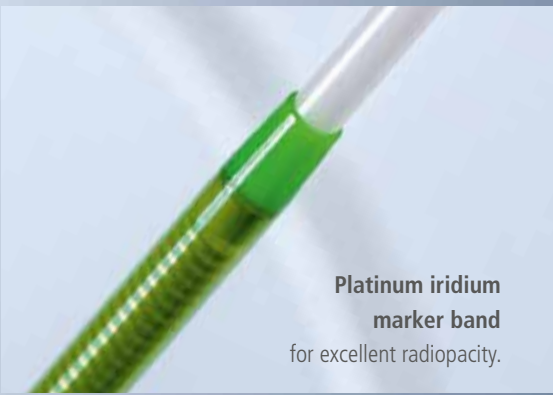

**Platinum iridium
marker band**
for excellent radiopacity.

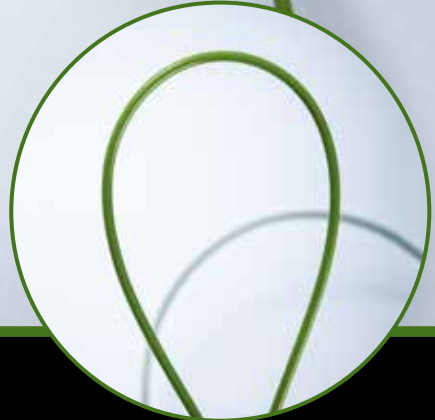

Designed for exceptional **🎯**
kink resistant performance
and non-ovaling effect.

2020 Contract Medical International. All rights reserved. DS-0004.05.2020

CMI 🎯 A product of
Contract Medical International GmbH
Lauensteiner Str. 37
01277 Dresden • Germany
Tel. +49 351 213 88 88
HeraeusMedicalComponents.com

Contract Medical International
joins Heraeus Medical Components
CMI 🎯 | Heraeus

DuraSheath
Kink resistant introducer sheath

🎯 Rethink the basics


Supporting complex vascular access

DuraSheath




Kink resistant introducer sheath

From easy insertion to smooth navigation to superb visibility: DuraSheath is designed to combine optimal performance with excellent pushability and trackability.

While maintaining the inner diameter, DuraSheath allows for smooth device passage with an exceptionally small kink radius and non-ovaling effect. The detachable hemostatic valve provides flexibility and versatility, allowing you the ability to focus on the intervention.

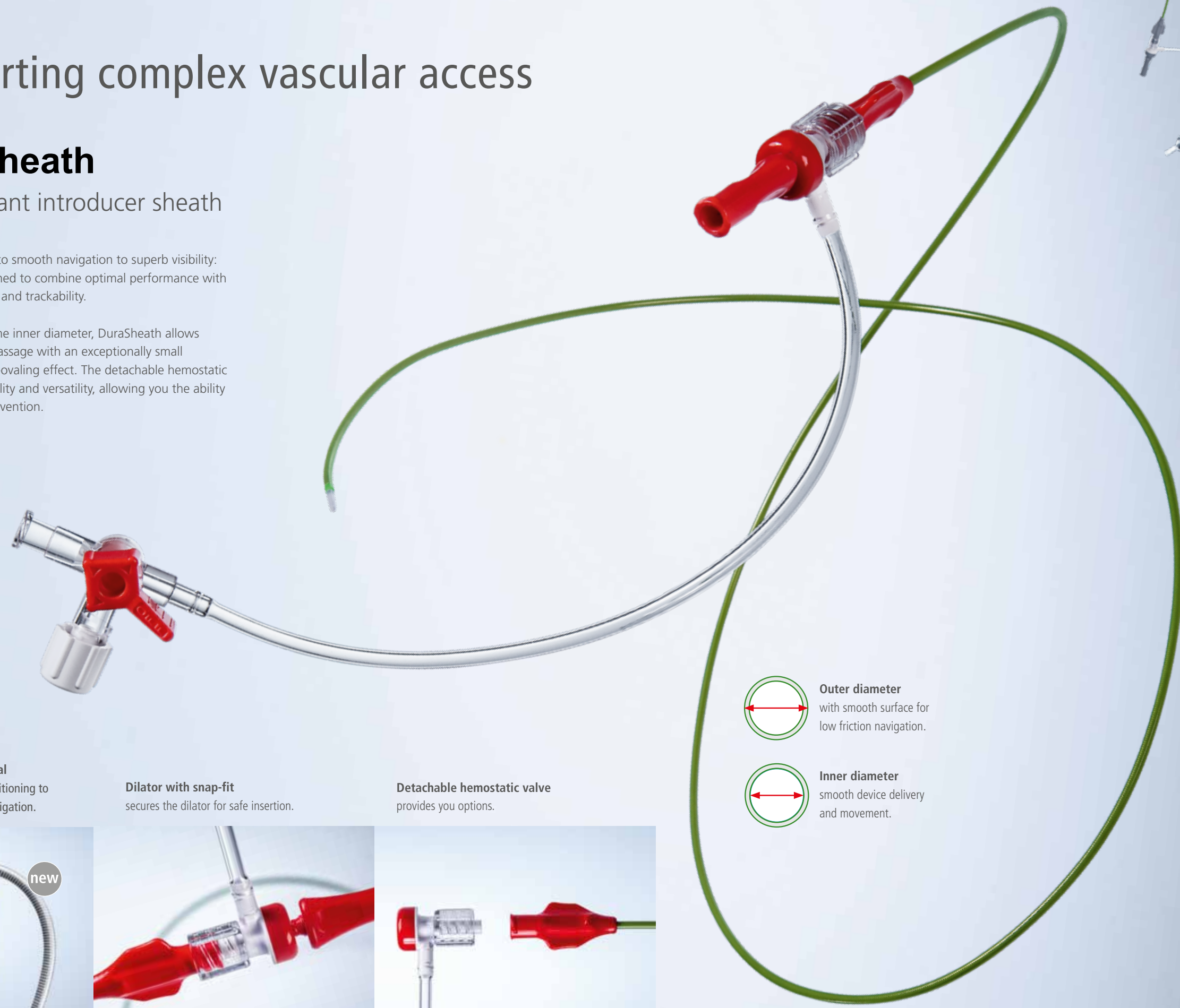

new Optional tip designs

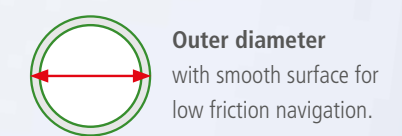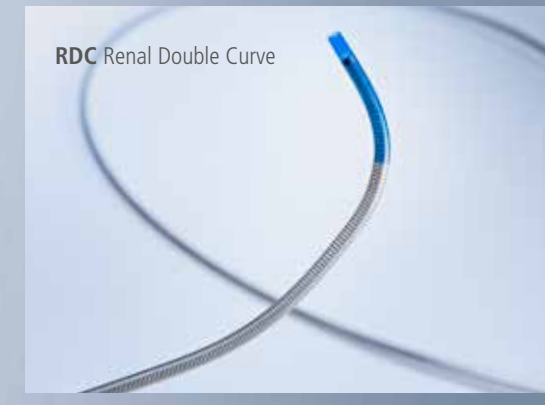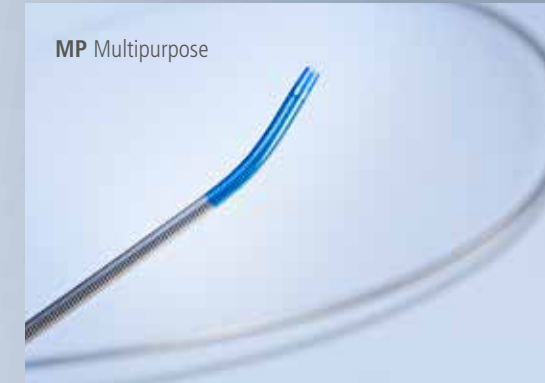

New optional shaft material
optimal backup support transitioning to softer distal zone for safe navigation.



Dilator with snap-fit
secures the dilator for safe insertion.



Detachable hemostatic valve
provides you options.




DuraSheath

Technical details / Ordering information

Providing access to target vessel not limited to

- peripheral
- renal
- coronary
- carotid

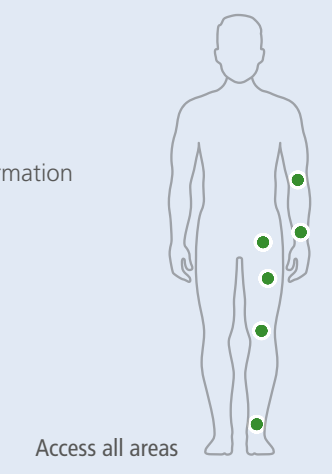

Introducer Sheath	
Material – outer	Stainless steel reinforced polymer
Material – inner	PTFE liner
Flexibility	Regular & Firm
Sheath coating	Hydrophobic
Sheath tip	Smooth taper with radiopaque platinum iridium marker band
French size	4, 5, 6, 7, 8 Fr
Length	15, 45, 60, 90, 130 cm
Sidearm/hemostatic valve	Detachable, with color coded 3-way-stopcock

Tip design	Tip design						Inner Diameter Min. (mm/inch)	Units/Case
	Straight	Straight	Multipurpose MP	Left Int. Mammary Artery LIMA	Renal Double Curve RDC	Hockey Stick HS		
Shaft version	Regular	Firm	Firm	Firm	Firm	Firm		
Diam. (Fr)	4 Fr	4 Fr	4 Fr	4 Fr	4 Fr	4 Fr	1.60 mm / 0.063"	5
Length (cm)	15 cm	45 cm	60 cm	90 cm	130 cm	15 cm	1.60 mm / 0.063"	5
Catalogue Number	FG-04015	FG-04045	FG-04060	FG-04090	FG-04130-F		1.60 mm / 0.063"	5
	5 Fr	5 Fr	5 Fr	5 Fr	5 Fr	5 Fr	1.92 mm / 0.075"	5
	15 cm	45 cm	60 cm	90 cm			1.92 mm / 0.075"	5
		FG-05045-F	FG-05045-F-MP	FG-05045-F-LIMA	FG-05045-F-RDC	FG-05045-F-HS	1.92 mm / 0.075"	5
	6 Fr	6 Fr	6 Fr	6 Fr	6 Fr	6 Fr	2.25 mm / 0.088"	5
	15 cm	45 cm	60 cm	90 cm	130 cm		2.25 mm / 0.088"	5
		FG-06045-F	FG-06045-F-MP	FG-06045-F-LIMA	FG-06045-F-RDC	FG-06045-F-HS	2.25 mm / 0.088"	5
	7 Fr	7 Fr	7 Fr	7 Fr	7 Fr	7 Fr	2.55 mm / 0.100"	5
	45 cm	60 cm	90 cm				2.55 mm / 0.100"	5
		FG-07045-F	FG-07045-F-MP	FG-07045-F-LIMA	FG-07045-F-RDC	FG-07045-F-HS	2.55 mm / 0.100"	5
	8 Fr	8 Fr	8 Fr	8 Fr	8 Fr	8 Fr	2.88 mm / 0.113"	5
	45 cm	60 cm	90 cm				2.88 mm / 0.113"	5
		FG-08045-F	FG-08045-F				2.88 mm / 0.113"	5
	8 Fr	8 Fr	8 Fr				2.88 mm / 0.113"	5
	60 cm	90 cm					2.88 mm / 0.113"	5
		FG-08060-F	FG-08060-F				2.88 mm / 0.113"	5
	8 Fr	8 Fr	8 Fr				2.88 mm / 0.113"	5
	90 cm						2.88 mm / 0.113"	5
		FG-08090-F	FG-08090-F				2.88 mm / 0.113"	5

15 cm versions contain two dilators.

Dilator Guidewire Compatibility	4 Fr / 5Fr / 6Fr		7 Fr / 8 Fr		
	Length	15 cm	45 cm / 60 cm / 90 cm / 130 cm	15 cm	45 cm / 60 cm / 90 cm / 130 cm
Dilator Guidewire Comp. max.		0.035" & 0.018"	0.035"	0.038" & 0.018"	0.038"